


ARS6PV

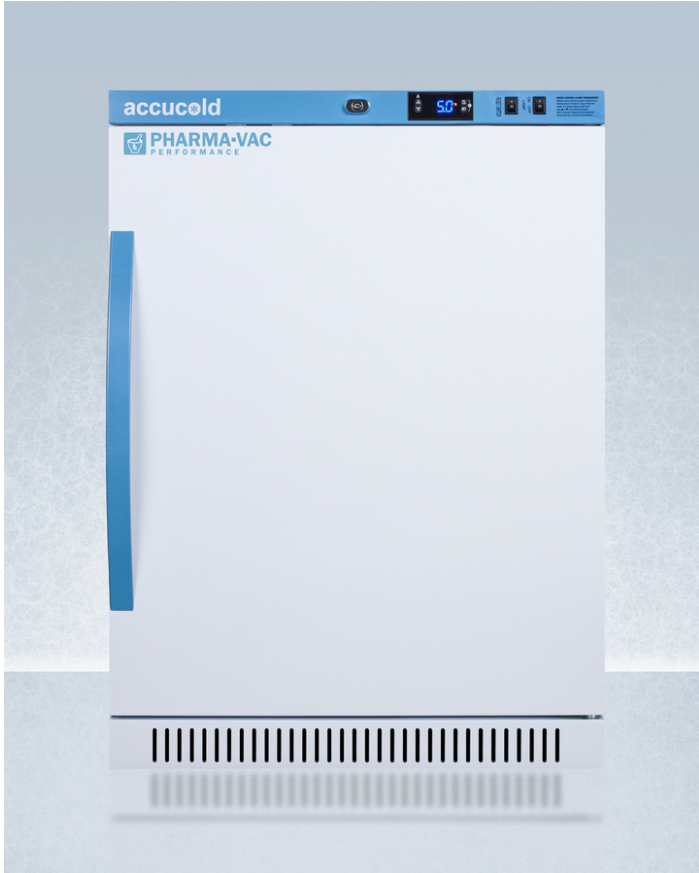

32.5" x 23.38" x 24.5" (H x W x D)

Performance series Pharma-Vac 6 cu.ft. freestanding ADA height all-refrigerator for vaccine storage

Highlights:

Intelligent microprocessor digital temperature controller

Improved temperature accuracy with a separate temperature probe for continuous product readout to simulate stored product temperature

Additional sensor reads interior air temperature for superior temperature control

Product Features:

Purpose-built design	Designed and purpose-built for pharmacy, medication, and vaccine applications to support meeting CDC/VFC guidelines
+/-1°C Variation	The published temperature variation is derived from the maximum deviation of an NTC sensor in a 1 oz. vial located nearest the chamber geometric center during a 24 hour test period
Adjustable temperature	Intuitive controls let you adjust the interior temperature between +2 to +8°C
Temperature alarm	Visual and audible alarms alert if the interior temperature rises or falls out of set range
Remote alarm contacts	Electronic controller is equipped with a volt-free contact for alarm forwarding to an external remote system (contacts located in the external rear of the unit)
Automatic defrost with adjustable cycle defrost	Adjustable cycle defrost operation ensures the unit avoid temperature spikes and ice buildup
Access port	Includes a pre-installed 1/2" probe hole in the right interior wall to allow easier connection to user-provided monitoring equipment
Open door alarm	Alarm will sound if the door is left ajar for over 4 minutes (adjustable from 1-99 minutes)
Audible power failure alarm	Battery-powered alarm sounds if the power connection is interrupted

ARS6PV Specifications:

Overview	
Height of Cabinet	32.5" (83 cm)
Width	23.38" (59 cm)
Depth	24.5" (62 cm)
Depth with Handle	25.5" (65 cm)
Capacity	6.0 cu.ft. (170 L)
Defrost Type	Automatic
Door	White
Cabinet	White
US Electrical Safety	ETL
Canadian Electrical Safety	ETL-C
Amps	0.75
Voltage/Frequency	115 V AC/60 Hz
Weight	120.0 lbs. (54 kg)
Shipping Weight	130.0 lbs. (59 kg)
Estimated Time to Ship	1 to 2 days
Parts & Labor Warranty	2 Years
Compressor Warranty	5 Years
Refrigerator Features	
Door Swing	RHD
Reversible	Yes
Shelf Type	Wire
Shelf Qty	4
Adjustable Shelves	Yes
Thermostat Type	Digital
Fan Type	Interior
Refrigerant Type	R600a
Refrigerant Amount	1.66oz.
Interior Light	Yes
Temperature Range	2 to 8°C

