

# PTHC125G



60" x 23.63" x 23.75" (H x W x D)

Single chamber passive heating blanket warmer with a self-closing glass door, stainless steel interior, digital thermostat, stainless steel cabinet, and lock

### Highlights:

Designed and assembled in the USA

No internal fans or moving parts ensure silent, safe, and efficient operation



Stainless steel construction inside and out for durable use and easy cleaning

### Product Features:

|                                 |  |
|---------------------------------|--|
| 90 to 140°F temperature range   | Designed and set to keep blankets warm in hospitals, medical and physical therapy offices, spas, and other facilities                  |
| Passive heating                 | Fixed heaters located at the bottom of the cabinet with no fans or moving parts to minimize lint buildup for safer, more efficient use |
| Digital thermostat              | Electronic controls are conveniently located for easy temperature adjustment   |
| Electronic readout              | Large LED temperature display ensures your staff can easily monitor the interior temperature   |
| Silent operation                | No moving parts to ensure a quiet performance  |
| Double pane tempered glass door | Heat-safe glass door offers a full view of stored contents, with a stainless steel door trim to highlight the interior                 |
| Professional handle             | Stainless steel handle completes the professional look   |
| Self-closing door               | Door will automatically close when left ajar, preventing items from accidentally being left exposed to ambient temperatures            |
| Stainless steel interior        | Interior walls are constructed from stainless steel to ensure lasting durability and easy cleaning                                     |

## PTHC125G Specifications:

| Overview                   |                      |
|----------------------------|----------------------|
| Height of Cabinet          | 60.0" (152 cm)       |
| Width                      | 23.63" (60 cm)       |
| Depth                      | 23.75" (60 cm)       |
| Depth with Handle          | 25.63" (65 cm)       |
| Depth with door at 90°     | 46.0" (117 cm)       |
| Capacity                   | 10.57 cu.ft. (299 L) |
| Door                       | Glass                |
| Cabinet                    | Stainless Steel      |
| US Electrical Safety       | ETL                  |
| Canadian Electrical Safety | ETL-C                |
| Amps                       | 5.5                  |
| Voltage/Frequency          | 115 V AC/60 Hz       |
| Weight                     | 160.0 lbs. (73 kg)   |
| Shipping Weight            | 175.0 lbs. (79 kg)   |
| Parts & Labor Warranty     | 1 Year               |
| UPC                        | 761101066127         |
| Warming Cabinet            |                      |
| Door Swing                 | RHD                  |
| Reversible                 | No                   |
| Adjustable Shelves         | Yes                  |
| Shelf Type                 | Chrome               |
| Shelf Qty                  | 4                    |
| Thermostat Type            | Digital              |
| Interior Height            | 45.25" (115 cm)      |
| Interior Width             | 20.0" (51 cm)        |
| Interior Depth             | 16.5" (42 cm)        |
| Dolly                      | Optional             |
| Temperature Range          | 90 to 140°F          |

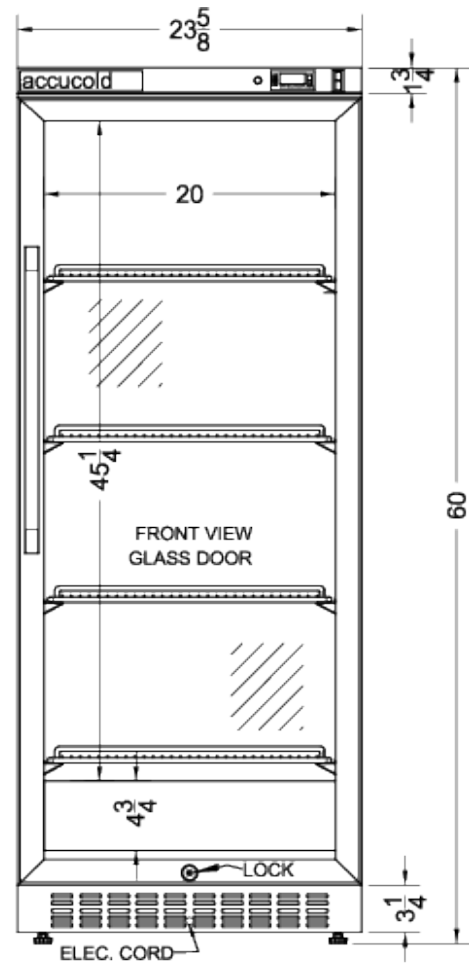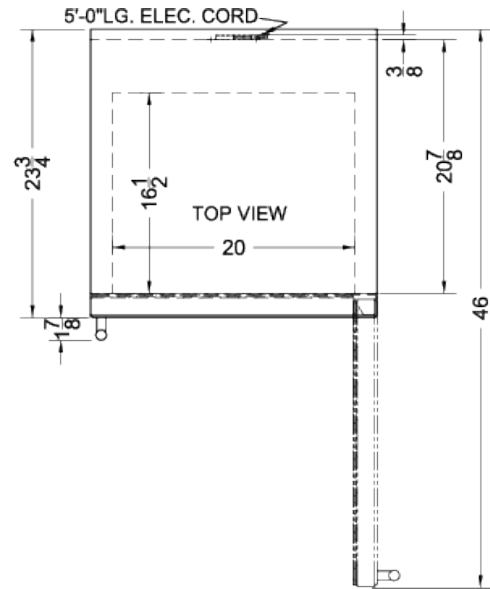

Puretherm Warming Cabinet