


ARS6MLDR

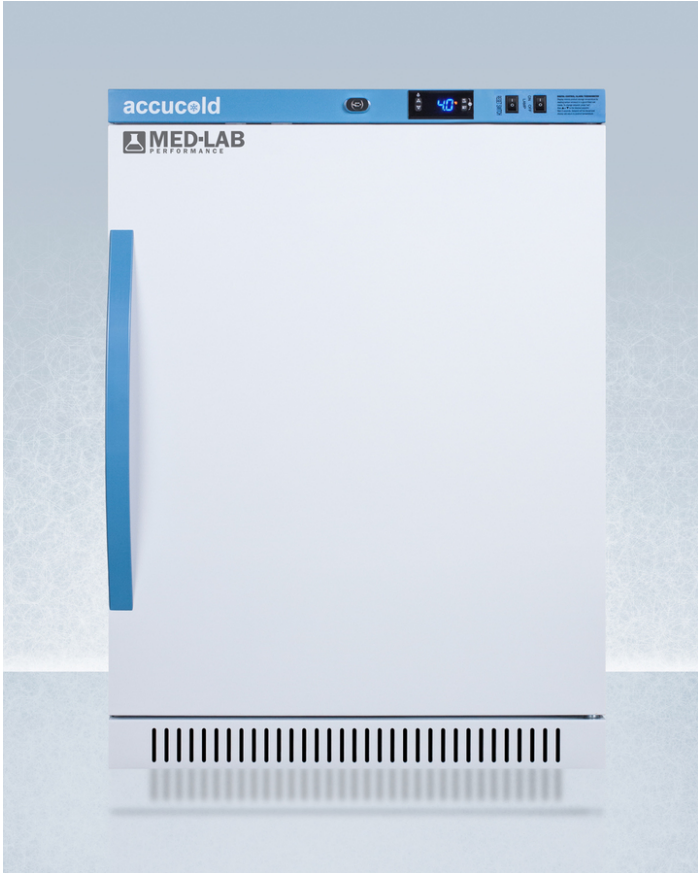

32.5" x 23.38" x 24.63" (H x W x D)

Performance series Med-Lab freestanding 6 cu.ft. ADA height all-refrigerator with three ventilated removable drawers, antimicrobial silver-ion handle, and hospital grade cord with 'green dot' plug

Highlights:

Powder-coated handle with naturally occurring silver ions to help reduce germ spread

Hospital grade cord with 'green dot' plug for added safety in high traffic areas

Intelligent microprocessor digital temperature controller

Product Features:

Purpose-built design	Designed and purpose-built to meet the demanding needs of scientific, life science, laboratory, and medical cold storage applications
+/-1°C Variation	The published temperature variation is derived from the maximum deviation of an NTC sensor in a 1 oz. vial located nearest the chamber geometric center during a 24 hour test period
Antimicrobial silver ion-coated handle	Powder-coated handle with naturally occurring silver ions to help reduce the spread of germs
Adjustable temperature	Intuitive controls let you adjust the interior temperature between +2 to +10°C
Full-extension ventilated drawers	Three drawers made from powder-coated aluminum glide out smoothly for convenient storage, with perforated bottoms to ensure proper air circulation
Removable drawers	Drawers can be easily removed for convenient transporting of contents to a lab counter or work station
Temperature alarm	Visual and audible alarms alert if the interior temperature rises or falls out of set range
Remote alarm contacts	Electronic controller is equipped with a volt-free contact for alarm forwarding to an external remote system (contacts located in the external rear of the unit)
Automatic defrost with adjustable cycle defrost	Adjustable cycle defrost operation ensures the unit avoid temperature spikes and ice buildup

ARS6MLDR Specifications:

Overview	
Height of Cabinet	32.5" (83 cm)
Width	23.38" (59 cm)
Width with Door Open	23.88" (61 cm)
Depth	24.63" (63 cm)
Depth with Handle	26.13" (66 cm)
Depth with door at 90°	46.5" (118 cm)
Capacity	6.0 cu.ft. (170 L)
Defrost Type	Automatic
Door	White
Cabinet	White
US Electrical Safety	ETL
Canadian Electrical Safety	ETL-C
Amps	0.75
Voltage/Frequency	115 V AC/60 Hz
Weight	130.0 lbs. (59 kg)
Shipping Weight	141.0 lbs. (64 kg)
Parts & Labor Warranty	2 Years
Compressor Warranty	5 Years
UPC	761101075129
Refrigerator Features	
Door Swing	RHD
Reversible	Yes
Shelf Type	Drawers
Shelf Qty	3
Adjustable Shelves	Yes
Thermostat Type	Digital
Fan Type	Interior
Refrigerant Type	R600a
Refrigerant Amount	1.66oz.
Interior Light	Yes
Temperature Range	2 to 10°C
Sound Level	36.4 dB
Dimensions	
Interior Height	24.5" (62 cm)
Interior Width	20.0" (51 cm)
Interior Depth	17.0" (43 cm)
Compressor Step Height	6.5" (17 cm)

Compressor Step Width	20.0" (51 cm)
Compressor Step Depth	5.0" (13 cm)

