


ARS6PV-AFZ5PVBIADASTACK

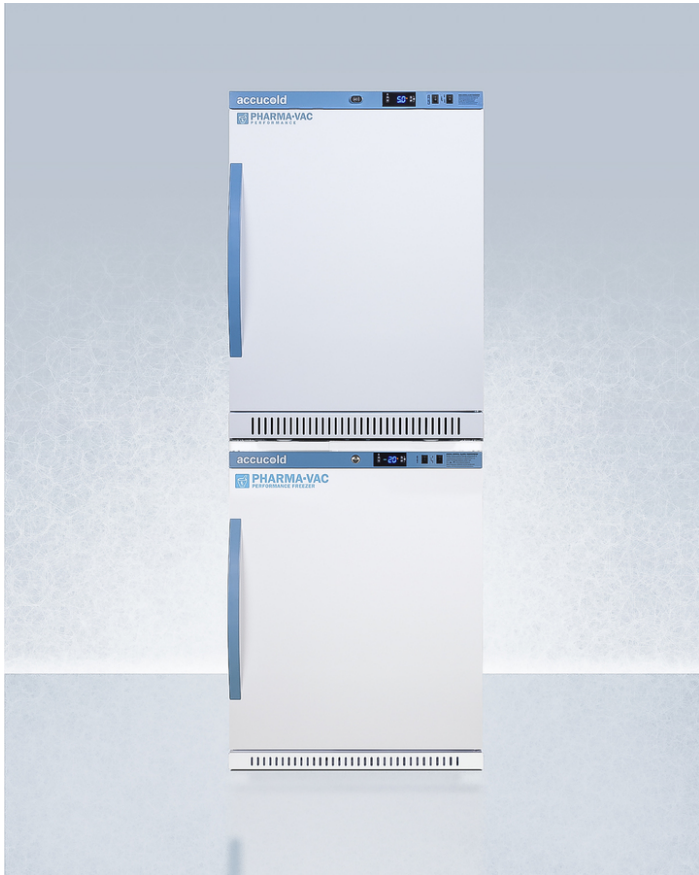

64" x 23.38" x 25" (H x W x D)

PHARMA-VAC Performance Series stacked combination of ARS6PV vaccine all-refrigerator and AFZ5PVBIADA manual defrost all-freezer

Highlights:

Purpose-built design to store medication, pharmaceuticals, and other sensitive materials

Independent refrigerator and freezer storage in one slim-fitting footprint for space-challenged facilities

Intelligent microprocessor digital temperature controller for optimum temperature management

Product Features:

Performance Series Refrigerator-Freezer Combination	Purpose-built all-refrigerator and all-freezer for vaccine storage in medical, pharmaceutical, and life science applications
Purpose-built design	Designed and purpose-built for pharmacy, medication, and vaccine applications to support meeting CDC/VFC guidelines
+/-1°C Variation	The published temperature variation is derived from the maximum deviation of an NTC sensor in a 1 oz. vial located nearest the chamber geometric center during a 24 hour test period
Space-saving design	Allows you to create a full refrigerator-freezer with independent controls in a slim-fitting footprint
Easy assembly	Unit ships in three boxes for customer to install the StackRack on the bottom refrigerator/freezer and place and secure the other unit on top
Antimicrobial silver ion-coated handles	Powder-coated handles with naturally occurring silver ions to help reduce the spread of germs
Adjustable temperature	Intuitive controls let you adjust the interior temperature between +2 to +8°C (refrigerator) and -30° to -10°C (freezer)
Temperature alarms	Visual and audible alarms alert if the interior temperature rises or falls out of set range

ARS6PV-AFZ5PVBIASTACK Specifications:

Overview	
Height of Cabinet	64.0" (163 cm)
Width	23.38" (59 cm)
Width with Door Open	23.88" (61 cm)
Depth	25.0" (64 cm)
Depth with Handle	26.75" (68 cm)
Depth with door at 90°	47.0" (119 cm)
Capacity	9.88 cu.ft. (280 L)
Door	White
Cabinet	White
US Electrical Safety	ETL
Canadian Electrical Safety	ETL-C
Voltage/Frequency	115 V AC/60 Hz
Shipping Weight	290.0 lbs. (132 kg)
Parts & Labor Warranty	2 Years
Compressor Warranty	5 Years
UPC	761101104133
Refrigerator-Freezer	
Door Swing	RHD
Reversible	Factory Reversible
Adjustable Shelves	Yes
Interior Light	Yes
Refrigerator - Interior Height	24.5" (62 cm)
Refrigerator - Interior Width	20.0" (51 cm)
Refrigerator - Interior Depth	17.0" (43 cm)
Refrigerator - Capacity	6.0 cu.ft
Refrigerator - Shelf Type	Wire
Refrigerator - Shelf Qty	4
Refrigerator - Temperature Range	2 to 8°C
Refrigerator - Refrigerant Type	R600a
Refrigerator - Refrigerant Amount	1.66 oz.
Freezer - Interior Height	22.5" (57 cm)
Freezer - Interior Width	18.0" (46 cm)
Freezer - Interior Depth	18.0" (46 cm)
Freezer - Capacity	3.88 cu.ft

Freezer - Shelf Type	Wire
Freezer - Shelf Qty	4
Thermostat Type	Digital
Fan Type	Interior
Freezer - Temperature Range	-30 to -10°C
Freezer - Refrigerant Type	R290a
Freezer - Refrigerant Amount	1.73 oz.

